


SINGLE USE SUPPORT.

PIONEERING BIOPHARMA

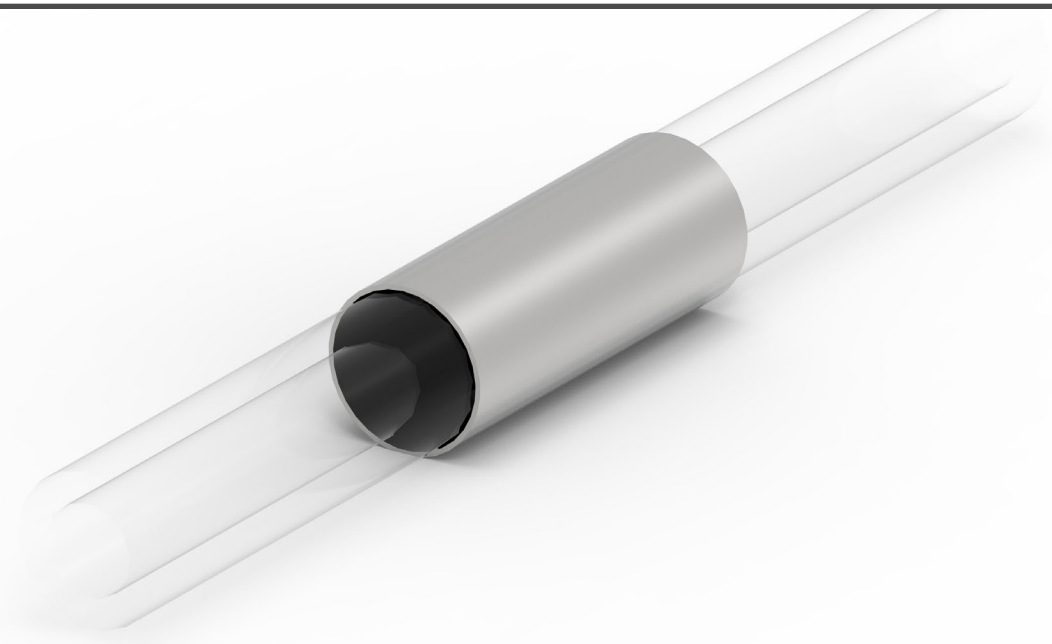

IRIS SEAL

Pinch pipe for crimping

*Sealing single-use assemblies.
Improved sterile process flexibility.
Compatible with all tubings.*


IRIS SEAL

Aseptic Disconnection.

IRIS Seals are metallic pinch pipes for crimping that can be pre-assembled onto single-use assemblies and allow easy aseptic disconnection while ensuring sterility and protection from the external environment.

Facts:

- IRIS seals are cylindrical metal pinch pipes allowing aseptic disconnection with crimping
- Withstanding pressure of up to 1.5 bar
- Available in standard outer diameter sizes of tubings: 1/2", 1/4" 5/8" - customization upon request
- Suitable for cleanroom and in compliance with international biopharmaceutical standards

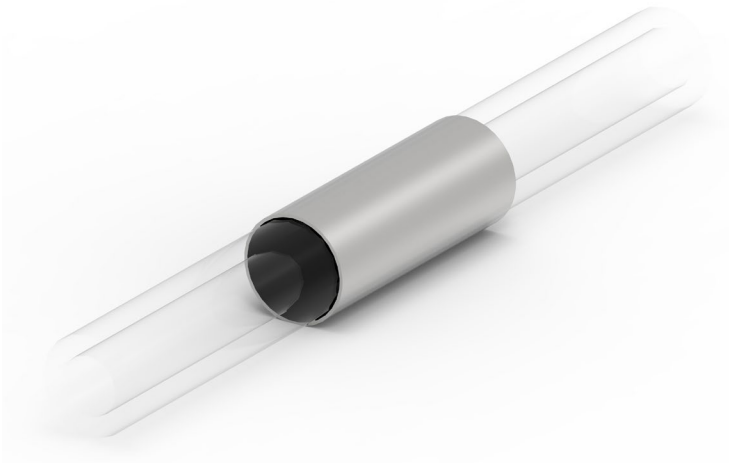


Single Use Support's Seals Crimping System

Aseptic Disconnection.

Single Use Support's Seals Crimping System provides a safe, reliable, and sterile disconnection of tubing in one single step. The light and powerful hand-held electrohydraulic tool allows to crimp and separate tubing with previously installed IRIS Seals and guarantees product integrity of the high-quality drug substances.

The system is compatible with various tubing sizes, pre-assembled with Single Use Support Seals, and requires one specific set of Single Use Support Dies per size, to seal and crimp the tube.

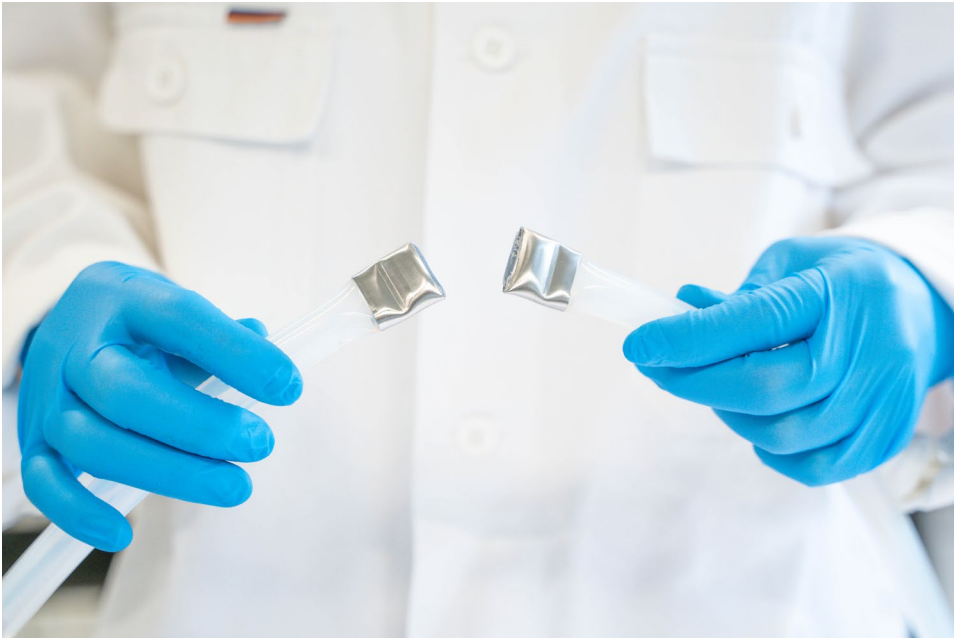

Single Use Support provides you with any highly qualified support all along the way. Contact us for further details!